


# Data Processing in MNova NMR

Below is a series of video tutorials covering different topics of MNova NMR processing. It is advisable to watch them in a sequence as they refer to previously discussed tools.

## Episode 1. Which MNova do I need?

This is the first episode of the series on NMR processing in MNova. Here I talk about MNova Lite CDE and full MNova NMR at Notre Dame campus.

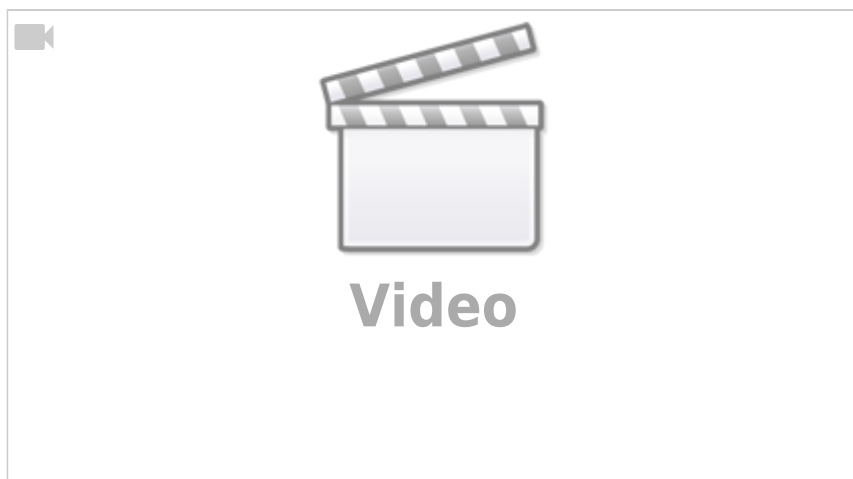


### Summary

- Use MNova Lite CDE for your 1D NMR work - unlimited licenses on campus
- Transition to MNova NMR (full version) for 2D NMR work - only a few licenses available
- Keep internet connection on while launching MNova and loading NMR data

---

## Episode 2. Creating MNova Document

This video covers basic toolbars and menus of the MNova window. We create a new document with a title page and save it.

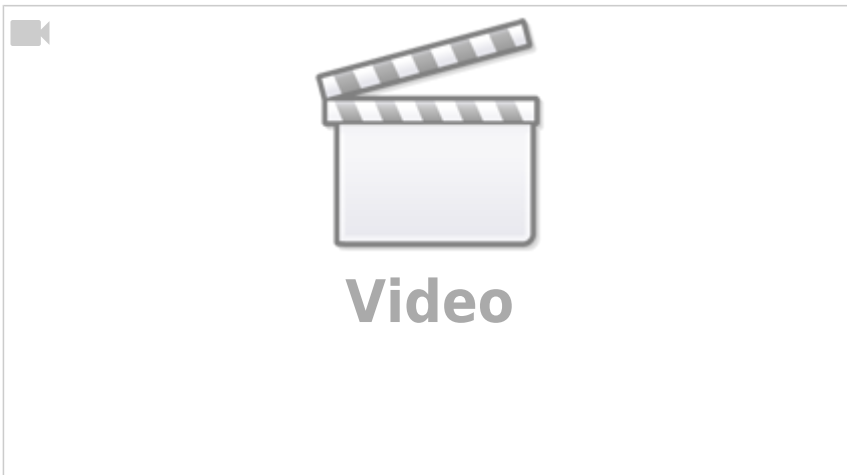


## Summary

- Make a title page with the basic sample info for each document
- Save with the original NMR data
- Save the MNova document frequently as you progress through your project

---

## Episode 3. Load data and move around with zoom tools

This episode demonstrates how to load NMR data into MNova document and explore the NMR spectrum using zoom tools.

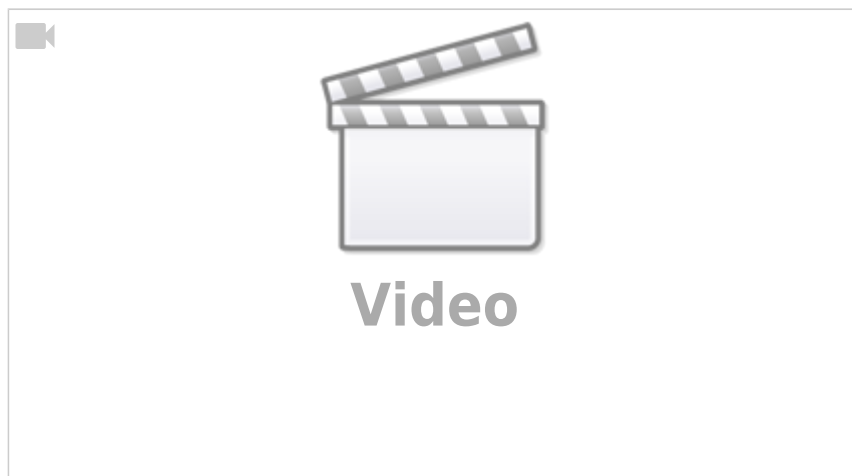

## Summary


- Load NMR data through drag-and-drop or Data Browser
- Update the spectrum Comment field and adjust its fonts size for display
- Add a title to the page itself to be able to find it easier in the Pages view
- Save frequently! (Cmd-S or Ctrl-S)
- Zoom intensity with a mouse wheel
- Zoom with 'Z', unzoom with Shift-'Z'; repeat a keystroke to switch the zooming direction
- Scale up with 'H' and 'F'. Pan with P.
- Navigate a spectrum with Shift-Z-Z-Z mode: click to zoom out, drag a box to zoom in

- Hit ESC on a keyboard to exit a zoom mode (or any other MNova tool)

---

## Episode 4. Manual Phase and Baseline Correction


This episode demonstrates a brief workflow for manual phase correction in a full version of MNova NMR.

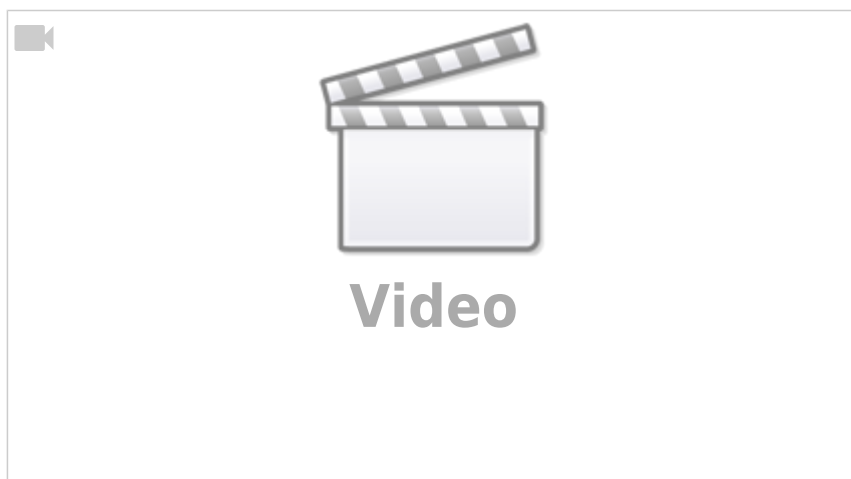

---

## 1D NMR Processing: Reference Deconvolution

This episode introduces a way to improve peak shapes in the spectrum if they were distorted by imperfect shimming. The tutorial focuses on basic concepts and a practical workflow. For more background on reference deconvolution, see this [blog post](#)

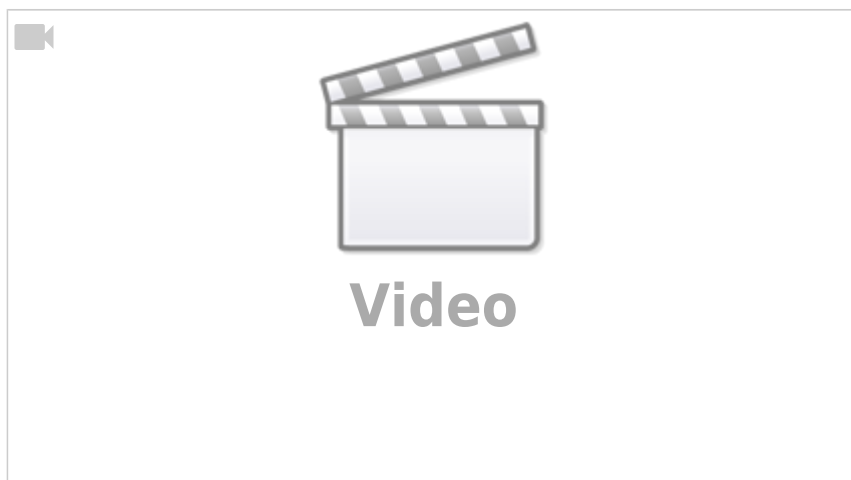

Work in progress! Future episodes are coming...

From:

<https://pydio.campus.nd.edu/docs/> - **NMR**

Permanent link:

[https://pydio.campus.nd.edu/docs/doku.php?id=mnova\\_data\\_processing&rev=1746372227](https://pydio.campus.nd.edu/docs/doku.php?id=mnova_data_processing&rev=1746372227)

Last update: **2025/05/04 16:23**

